



CE122

DIGITAL INTERFACE

Electronic interface that connects TecQuipment's Control Engineering (CE) range to a suitable computer for computer control and data acquisition



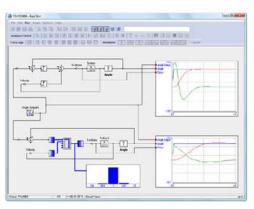

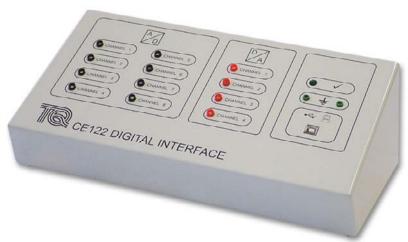

A TYPICAL SCREENSHOT OF THE CE2000 SOFTWARE (SUPPLIED)

- Compact, bench-mounting unit
- Easy to connect, low-voltage inputs and outputs are buffered for educational use
- Supplied with TecQuipment's Control Software (CE2000) for control and data acquisition
- Fully digital needs no user adjustments
- Can work as an interface to other compatible laboratory equipment for computer control and data acquisition
- Connects to most modern computers



TECQUIPMENT LTD, BONSALL STREET, LONG EATON, NOTTINGHAM NGIO 2AN, UK TECQUIPMENT.COM +44 115 972 2611 SALES@TECQUIPMENT.COM

PE/db 0417 Page 1 of 2

DIGITAL INTERFACE

DESCRIPTION

The Digital Interface is for use with most products in TecQuipment's Control Engineering range, but it will also work with any other compatible laboratory equipment.

It is an alternative to the CE120 Controller, when the user only needs the interface part of the CE120.

It converts analogue inputs from other equipment into digital signals for a computer. It also converts the digital signals from a computer into analogue signals to control other equipment.

It will allow the user to monitor up to eight input signals and control up to four output signals.

The interface includes a copy of TecQuipment's CE2000 Control Software (see separate datasheet for details).

AVAILABLE EXPERIMENT MODULES

Selected products in TecQuipment's Control Engineering range:

- Thermal Process Control (CE103)
- Coupled Tanks (CE105 and CE105 MV)
- Ball and Beam (CE106)
- Engine Speed Control (CE107)
- Coupled Drives (CE108)
- Ball and Hoop (CE109)
- Servo Trainer (CE110)

STANDARD FEATURES

- Supplied with comprehensive user guide
- Five-year warranty
- Made in accordance with the latest European Union directives

ESSENTIAL ANCILLARIES

Suitable computer with two spare USB connections (not supplied by TecQuipment) – refer to CE2000 datasheet for full details.

ESSENTIAL SERVICES

ELECTRICAL SUPPLY:

90 VAC to 250 VAC 65 mA, 50/60 Hz, with earth

BENCH SPACE NEEDED:

400 mm x 300 mm

OPERATING CONDITIONS

OPERATING ENVIRONMENT:

Laboratory

STORAGE TEMPERATURE RANGE:

-25°C to +55°C (when packed for transport)

OPERATING TEMPERATURE RANGE:

+5°C to +40°C

OPERATING RELATIVE HUMIDITY RANGE:

80% at temperatures < 31°C decreasing linearly to 50% at 40°C

SOUND LEVELS

Less than 70 dB(A)

SPECIFICATIONS

NETT DIMENSIONS AND WEIGHT:

360 x 170 x 125 mm, 3 kg

PACKED DIMENSIONS AND WEIGHT:

0.04 m³, 4.7 kg (approx – packed for export)

INPUT AND OUTPUTS:

Eight A-D inputs: ±10 V Four D-A outputs: ±10 V

Type B USB socket connection for a suitable computer

OTHER PARTS INCLUDED:

Connecting cables



Page 2 of 2

PE/db 0417